medela 
mother's milk,
everyday amazing™

Sonata®

EN Instructions for use

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Read all instructions before using this product.
Review the most current product instructions at
Lea todas las instrucciones antes de usar este producto.
Lea las instrucciones más actualizadas del producto en
Veuillez lire toutes les directives avant d'utiliser ce produit.
Consultez les directives les plus récentes sur le produit à l'adresse
www.MedelaBreastfeedingUS.com | www.medela.ca

*All parts that touch breast milk, are made without BPA.
*Todas las partes que tocan la leche materna están fabricadas sin BPA.
*Toutes les pièces qui entrent en contact avec le lait maternel ne contiennent pas de BPA.

PLEASE SAVE THESE INSTRUCTIONS
CONSERVE ESTAS INSTRUCCIONES
VEUILLEZ CONSERVER CES DIRECTIVES

1. Important Safeguards

When using electrical products, especially when children are present, basic safety precautions should always be followed.

2

READ ALL INSTRUCTIONS BEFORE USING THIS PRODUCT.

The defined signal words identify all instructions that are important to safety. Failure to observe these instructions can lead to injury or damage to the breast pump or yourself. When used in conjunction with the following words, the defined signal words stand for:

WARNING Can lead to serious injury or death.

CAUTION Can lead to minor injury.

NOTE Can lead to material damage.

TIP Useful or important information that is not related to safety.

WARNING: To avoid fire, electrocution, or serious burns:

- This is a single-user product. Use by more than one person may present a health risk and voids the warranty.
- Always unplug electrical product immediately after use except for when charging.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- The breast pump and accessories are not heat-resistant: keep away from heated surfaces or open flames.
- Do not use near flammable materials.
- The breast pump should never be left unattended when plugged into a power source.
- Never operate an electrical device if it has a damaged cord or plug, is not working properly, or if it has been dropped or damaged. If damage is found, immediately discontinue use and call Medela Customer Service at 1-800-435-8316.
- Do not use an electrical outlet device that has been exposed to water or other liquids including:
 - Do not use while bathing or showering.
 - Never place or drop into water or other liquids.
 - Do not run water over breast pump.
 - If a device has been exposed to water or other liquids, do not touch, unplug the device from electrical outlet, turn off and contact manufacturer.

**WARNING: To avoid health risk and reduce the risk of injury:**

- This device cannot be serviced or repaired. Do not repair yourself. Do not modify the device.
- Never use a damaged device. Replace damaged or worn parts.
- Use only the power adaptor that comes with the breast pump.
- Use the breast pump only for its intended use as described in this manual.
- Do not use the breast pump while sleeping or overly drowsy.
- Do not use Sonata breast pump while operating a moving vehicle.
- Inspect all appropriate pump components before each use.
- If tubing becomes moldy, discontinue use and replace tubing.
- Pumping can induce labor. Do not pump until after giving birth. If you become pregnant while breastfeeding or breast pumping, consult with a licensed healthcare professional before continuing.
- If infected with Hepatitis B, Hepatitis C, or Human Immunodeficiency Virus (HIV), pumping breast milk will not reduce or remove the risk of transmitting the virus to your baby through your breast milk.
- Do not microwave or boil breast milk. Microwaving can cause severe burns to baby's mouth from hot spots that develop in the milk during microwaving. (Microwaving can also change the composition of breast milk.)
- Clean and sanitize all parts that come into contact with your breast and breast milk prior to first use.
- Wash all parts that come into contact with your breast and breast milk after every use.
- Only use Medela recommended parts with your Sonata® Breast Pump.
- Close supervision is necessary when the breast pump or accessories are used in the vicinity of children.

1. Important Safeguards (cont.)

4

CAUTION: Can lead to minor injury:

- Portable and mobile radio frequency communications equipment can affect the breast pump.
- Make sure the voltage of the power adaptor is compatible with the power source. See **Section 23** for technical specifications.
- Wash hands thoroughly with soap and water before touching breast pump, kit and breasts, and avoid touching the inside of bottles or lids.
- Separate and wash all parts that are exposed to breast milk immediately after use. This will help remove breast milk residue and prevent growth of bacteria.
- Always inspect breast shields, connectors, valves, membranes, bottles, lids, and tubing prior to use for cleanliness. Contact Medela Customer Service if cleaning does not resolve the issue.
- Only use drinking-quality tap or bottled water for cleaning your breast pump and parts.
- Do not store wet or damp parts as mold may develop.
- Do not run pump with wet tubing. Doing so may damage the breast pump.
- If you are experiencing discomfort at the base of the nipple due to rubbing of your breast tissue against the breast shield tunnel, use of a lubricant such as Tender Care™ Lanolin may be beneficial. For assistance with correct breast shield sizing and comfort please visit www.MedelaBreastshields.com or see a lactation consultant/breastfeeding specialist.
- If pumping is uncomfortable or causing pain, turn the unit off, break the seal between the breast and the breast shield with your finger and remove the breast shield from your breast.
- Contact your health care professional or breastfeeding specialist if you can express only minimal or no milk or if expression is painful. See **Section 7 & 8** for more information.
- While some discomfort may be felt when first using a breast pump, using a breast pump should not cause pain. For assistance with correct breast shield sizing and comfort please visit www.MedelaBreastshields.com or see a lactation consultant/breastfeeding specialist.
- Do not try to express with vacuum that is too high and uncomfortable (painful). The pain, along with potential breast and nipple trauma, may decrease milk output.
- Make sure tubing is not kinked or pinched while pumping.
- Do not hold the pump kit by the bottle. This can lead to blockage of the milk ducts and engorgement. See **Section 11** for more information.
- The breast pump uses Bluetooth® technology. Please follow your air carrier rules for the use of Portable Electronic Devices when flying with your breast pump.
- Using a breast pump on an aircraft is not recommended. The cabin pressure may affect the breast pumps performance.

Note

- Plastic bottles and parts become brittle when frozen and may break when dropped.
- Bottles and parts may become damaged if mishandled (e.g., dropped, over-tightened, or knocked over).
- Take appropriate care in handling bottles and components.
- Do not use the breast milk if bottles or components become damaged.

Indications for use

The Sonata® breast pump is a powered breast pump to be used by lactating women to express and collect milk from their breast. The Sonata Breast Pump is a single user device.

Contraindications for use

There are no known contraindications for use with this product.

Tip

If you have medical reasons or other needs for exclusively pumping, it is recommended that you use a hospital-grade (multi-user) breast pump such as the Medela Symphony® Breast Pump. To learn more, visit www.medelabreastfeedingus.com | www.medela.ca.

Product Description

This breast pump is a personal-use electric breast pump that includes 2-Phase Expression® technology and is capable of single and double pumping.

The operating life of this breast pump is defined to be approximately three 15-minute sessions per day, for one year. The operating life for the breast pump kit is 6 months.

**PLEASE SAVE THESE INSTRUCTIONS.
THIS IS A SINGLE-USER PRODUCT.**

Table of Contents

1. Important Safeguards	2-5
2. Product Description	8-9
3. Getting Started	10
4. Cleaning Overview	11
5. Cleaning Instructions	12-13
6. Breast Pump & Tubing Care	14-15
6A Cleaning the Breast Pump	14
6B Tubing Care	15
7. Breastfeeding Information	16
8. Breast Shield Sizing	17
9. Assembling Your Breast Pump	18-19
10. Powering Your Breast Pump	20-22
10A Charging Your Breast Pump	20
10B Battery Care & Status Codes	20-21
10C Power Adaptor Care	22
10D Traveling Outside of Your Country	22
11. Operating Your Breast Pump	23-26
11A Glossary of Terms, Controls & Icons	23
11B Single Pumping	24
11C Double Pumping	24
11D Pumping Steps	25-26
12. Finding Your Maximum Comfort Vacuum™	27
13. Features	28-33
Breast Pump Connectivity	28-29
Selecting a Rhythm	30
Setting the Session Timer	31
Pausing Your Breast Pump	32
Silencing Your Breast Pump Chimes	33

14. Storing Breast Milk.....	34
15. Preparing & Feeding Breast Milk.....	35
16. Troubleshooting.....	36-41
Air Leak Alert.....	36
Battery Drained.....	37
Battery Broken.....	38
Bluetooth® Broken.....	39
General Alert.....	40
Additional Troubleshooting.....	41
17. Long-term Storage.....	42
18. Disposal.....	43
19. Supplemental Information.....	44
20. Warranty.....	45
21. Meaning of Symbols.....	46
22. EMC/Technical Description.....	47-53
23. Technical Specifications.....	54

2. Product Description

8

Your breast pump system includes

**2x PersonalFit Flex™
Breast Shields
(21 mm)
Applied Part**

**2x PersonalFit
Breast Shields
(24 mm)
Applied Part**

For other breast shield sizes,
see **Section 8**

Replacement 101036653 (21 mm)
Replacement 101036654 (24 mm)

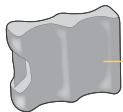
4x Bottles
Item 6007135

1x Bag
Item 3007396
Replacement 68052

1x Cooler bag
Item 3007304
Replacement 67068
(replacement includes ice
pack & 4 bottles with lids)


1x Ice pack
Item 8117010
Replacement 87092

Quality & performance
you can trust

**2x Breast shield
assemblies**

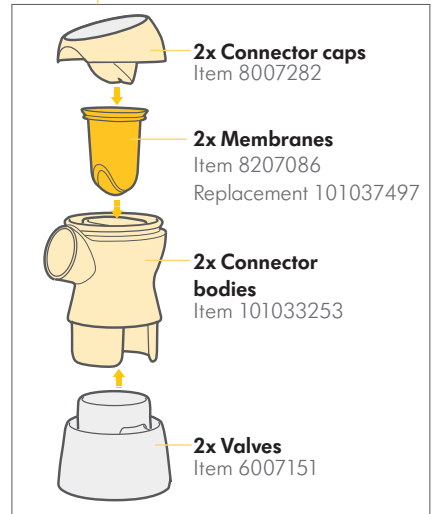Replacement 101037321


2x Connector caps
Item 8007282


2x Membranes
Item 8207086
Replacement 101037497

**2x Connector
bodies**
Item 101033253

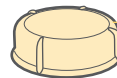
2x Valves
Item 6007151

4x Lids
Item 8107183

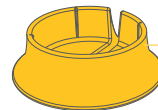
2x Bottle stands
Item 8107322

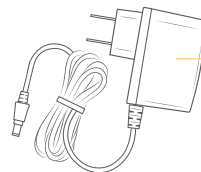(1x) Power adaptor
Item 101034947
Replacement
101037316
68050



Medela reserves the right to substitute any component or
accessory with a replacement of equivalent performance.

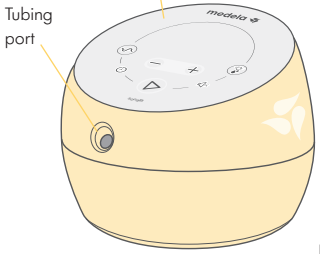1x Breast pump


1x Tubing
Item 8007278
Replacement
101037315

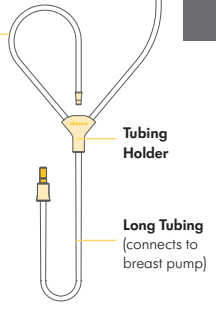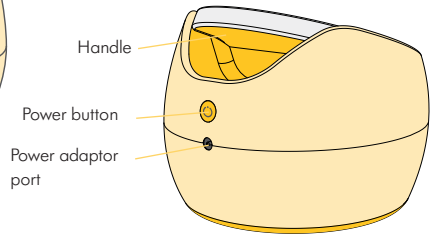
Short Tubing
(connects to breast shield assemblies)

Controls

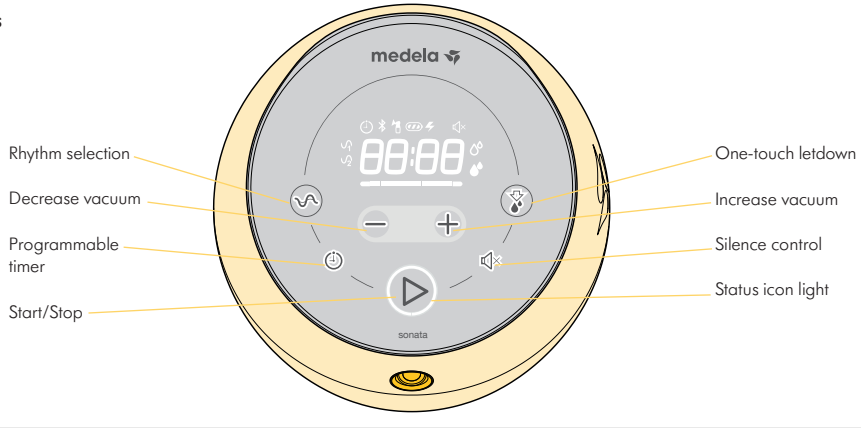Display Icons

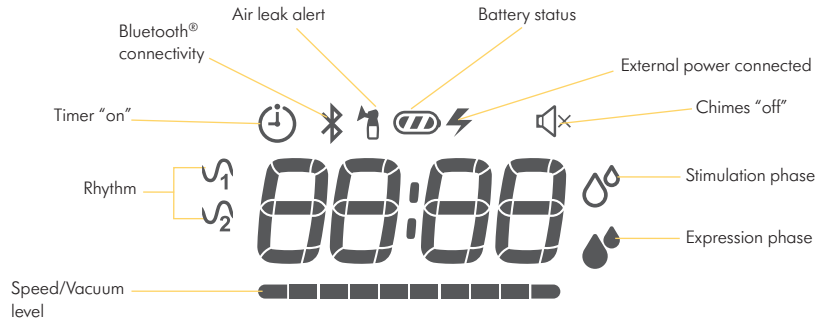3. Getting Started

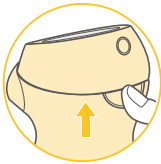
Note

It is important for you to do the following before using your breast pump for the first time:

1. Plug the Sonata into an external power source. Charge the battery for 12 hours. (See **Section 10**, Powering Your Pump).
2. Separate all parts that are exposed to the breast or breast milk prior to cleaning.
3. Clean - see **Section 5**.

Parts to clean

Separate all parts that are exposed to the breast and breast milk prior to cleaning.

Remove and dispose of the protective cover from the breast shield before first use.

Use your thumb to separate the connector cap and body.

Pull down on valve to remove from connector body. Do not rotate.

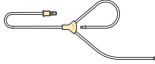
4. Cleaning Overview

Washing and sanitizing are two different activities. They must be done separately to protect you, your baby and the performance of your breast pump.

Wash – To clean the surfaces of the parts by physically removing contamination.

Sanitize – To kill living organisms, such as bacteria or viruses, that may be present on the surfaces of the parts.

When to Wash	Breast Shields 	Breast Pump Kit 	Bottles & Lids 
Before 1st use	<ul style="list-style-type: none"> ✓ wash ✓ sanitize 	<ul style="list-style-type: none"> ✓ wash ✓ sanitize 	<ul style="list-style-type: none"> ✓ wash ✓ sanitize
After each use	<ul style="list-style-type: none"> ✓ wash 	<ul style="list-style-type: none"> ✓ wash 	<ul style="list-style-type: none"> ✓ wash
Once per day	<ul style="list-style-type: none"> ✓ sanitize 	<ul style="list-style-type: none"> ✓ sanitize 	<ul style="list-style-type: none"> ✓ sanitize

When to Wash	Tubing 	Breast Pump 	Carry Bag 
As needed	<ul style="list-style-type: none"> ✓ See tubing care - Section 6 	<ul style="list-style-type: none"> ✓ wipe with clean, damp cloth 	<ul style="list-style-type: none"> ✓ wipe with clean, damp cloth

Detailed cleaning instructions on pages 12-15:

- Wash – Section 5
- Sanitize – Section 5
- Breast Pump and Tubing Care – Section 6

5. Cleaning Instructions

Supplies needed:

- Mild dish soap or Medela Quick Clean™ Breastmilk Removal Soap
- Clean dish towel or soft brush
- Clean wash basin
- Drinking-quality water
- Tongs
- Clean pot for boiling water

Parts to wash or sanitize:

- Breast shields
- Breast milk bottles
- Lids
- Connector bodies
- Valves
- Membranes
- Connector caps

Caution

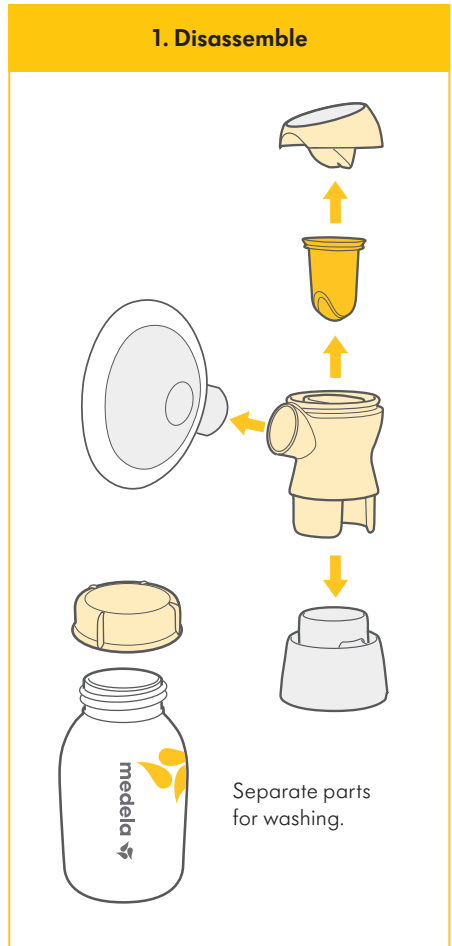
- Separate and wash all parts that are exposed to breast milk immediately after use. This will help to remove breast milk residue and prevent growth of bacteria.
- Only use drinking-quality tap or bottled water for cleaning.
- When cleaning the valves do not use small objects such as brushes. Be careful not to puncture.

Note

- Take care not to damage parts of the breast pump kit during cleaning.
- Store the dry breast pump kit in a clean bag/container until next use.
- Wash hands thoroughly.
- If you notice a white residue on your breast pump parts after boiling, your water may have a high mineral content. Remove residue by wiping parts with a clean cloth and allow to air dry.
- Distilled water is recommended when boiling parts to prevent substantial mineral build-up over time, which may compromise your parts.

1. Disassemble


Tip

- If using the dishwasher, parts may become discolored. This will not impact part function.
 - Tubing should only be washed if dirty or milk is present.
- See **Section 6** for more details.

2. Wash

Before first use and after each use

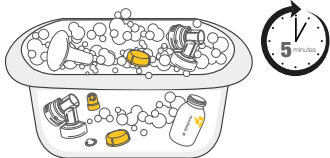
Sink/clean wash basin

- Step 1** – Rinse all separated parts that come into contact with the breast and breast milk in cool water to remove breast milk residue. Once rinsed, parts that have touched the breast or breast milk should be placed in a clean wash basin, dedicated to only cleaning these items. Do not place parts directly in the sink.
- Step 2** – Soak in warm soapy water for 5 minutes. Wash with a clean, unused dish cloth and rinse with clean water.
- Step 3** – Place parts on a clean surface and/or towel. Allow parts to air dry. Store dry parts in a cool place when not in use. Do not store wet or damp parts.

OR

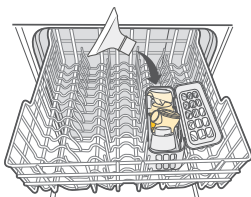Dishwasher

- Step 1** – Wash all separated parts by using the top rack of your dishwasher.
- Step 2** – Place parts on a clean surface and/or towel. Allow parts to air dry. Store dry parts in a cool place when not in use. Do not store wet or damp parts.

3. Sanitize

Before first use and once per day

Stovetop

- Step 1** – Wash hands thoroughly.
- Step 2** – Fill pot with enough water to cover all parts. Bring to a boil.
- Step 3** – Place parts in boiling water for 10 minutes.
- Step 4** – Allow parts to air dry. Store dry parts in a cool place when not in use. Do not store wet or damp parts.

OR

Microwave

- Step 1** – Use Medela's Quick Clean™ Micro-Steam™ bags. Sold separately (Follow instructions as provided on bag)
- Step 2** – Place parts on a clean surface and/or towel. Allow parts to air dry. Store dry parts in a cool place when not in use. Do not store wet or damp parts.

6. Breast Pump & Tubing Care

Supplies needed:

- Mild dish soap or Medela Quick Clean™ Breast Milk Removal Soap
- Clean dish towel
- Clean sink or bowl
- Drinking-quality water

Parts needed:

- Tubing
- Breast pump

Tip

- Tubing should only be washed if dirty or milk is present. Follow instructions in **Section 6B**.
- There is no need to clean tubing if condensation inside the tube is from previous washings or atmospheric conditions. Follow drying instructions in **Section 6B**.

6A. Cleaning the Breast Pump

1

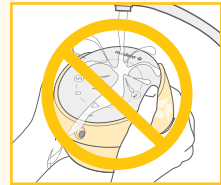
Wipe pump unit with a clean, damp cloth.

Warning

- Do not immerse the breast pump in water; do not run water over the breast pump.
- When cleaning the breast pump display, only use a damp, soft cloth with mild dish soap. Other cleaners or abrasive cloths may degrade the display performance.

6B. Tubing Care

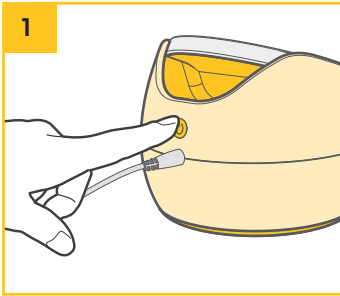
- Tubing does not need to be washed before first use.
- Inspect tubing after each pumping session.

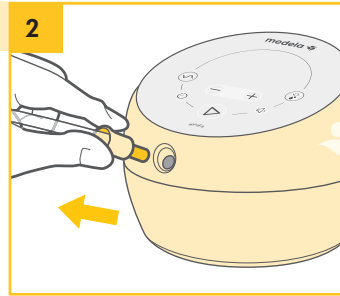
Caution

- Do not store wet or damp parts as mold may develop.
- If tubing becomes moldy, discontinue use and replace tubing. To find replacement parts visit www.ShopMedela.com (USA) or www.medelastore.ca (Canada)
- **Do not use breast pump with wet tubing; doing so will cause damage to your breast pump.**

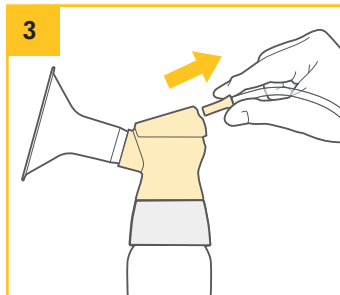
To properly clean your tubing


Turn off  breast pump and unplug from power source.

To remove tubing from the breast pump, grab the tube adaptor and pull directly from the tubing port. Do not wiggle or pull tubing at an angle.

Remove tubing from breast shield.

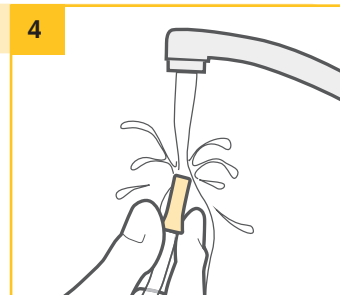- Rinse tubing in cool water.
- Wash tubing in warm, soapy water.
- Rinse tubing with clear water.
- Shake out water droplets.
- Hang to air dry.

7. Breastfeeding Information

Common pumping questions & answers

How often should you breast pump?

A breast pump is a replacement for when you are separated from your baby. It is important to pump when the baby would normally be breastfeeding. For example, a working mother may pump 2-3 times during an 8-hour work day.

How long should your pumping session last?

Pumping times can vary from mother to mother. Pumping sessions usually last from 15 to 30 minutes.

How should your breasts feel before and after pumping?

Before pumping, your breasts will have a firm, heavy feeling. After pumping, your breasts should feel soft with no firm areas. Firm areas could indicate that the breast is not completely drained.

8. Breast Shield Sizing

1. Medela's PersonalFit™ and PersonalFit Flex™ sizing

Pumping should not hurt. For maximum comfort and pumping efficiency, Medela offers PersonalFit Flex™ in four breast shield sizes*:

21 mm

24 mm

27 mm

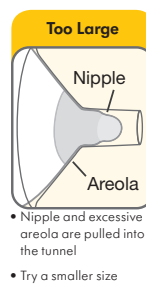
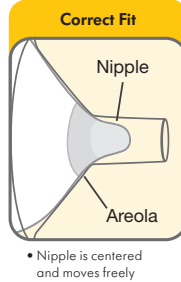
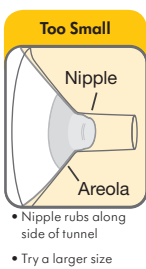
30 mm

* Medela's original PersonalFit breast shields are also available in 36 mm.

Visit MedelaBreastshields.com to view the full guide to help you determine your optimal size based on your nipple diameter.

2. Test your breast shield size

1. Start with the 21 mm that came with your pump, or the size determined by measuring.
2. Center nipple and gently hold breast shield against your breast.
3. Adjust for Maximum Comfort Vacuum™, to achieve optimum suction level.
4. Refer to images while pumping in expression.

3. Should you try a new size?

- Does your nipple rub sides of tunnel, to the point of causing discomfort?
- Do you see excessive areola being pulled into tunnel?
- Do you see any redness?
- Is your nipple or areola turning white?
- Do you feel unexpressed milk after pumping?

If you answered "YES" to any of these questions, consider trying a new size by following the measuring instructions above.

If you are still unsure if you selected the correct size, see a lactation consultant, breastfeeding specialist, or visit MedelaBreastshields.com for assistance in choosing the right size breast shield.

Caution

- While some discomfort may be felt when first using a breast pump, using a breast pump should not cause pain. If you are unsure about breast shield sizing, please contact a health care professional or breastfeeding specialist who can help you get a proper fit.
- If you are experiencing discomfort at the base of the nipple due to rubbing of your breast tissue against the breast shield tunnel, use of a lubricant such as Tender Care™ Lanolin may be beneficial.

Questions? Visit www.medela.com or call us at 1-800-435-8316.

9. Assembling Your Breast Pump

Parts Needed:

- Tubing
- Breast pump
- Breast shields
- Breast milk bottles
- Lids
- Connector bodies
- Valves
- Membranes
- Connector caps

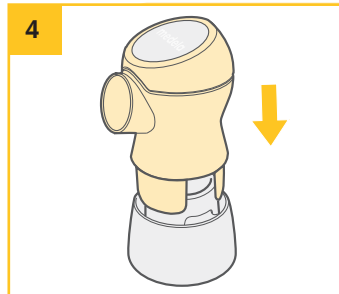
T Tip

- Only use Authentic Medela Spare Parts. See **Section 2** for details.
- Check breast pump kit components for wear or damage before use. Replace if necessary.
- Always inspect all parts prior to use for cleanliness. If dirty, see **Section 4**.
- To prevent damage to the breast pump all components must be completely dry before use.


T Tip

Wash hands thoroughly with soap and water before touching breast pump, breast pump kit and breasts. Avoid touching the inside of containers and lids.

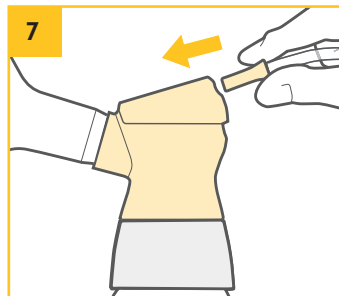
Using the two legs as guides, insert the connector body onto the valve.


Insert short tubing into the connector cap until snug.

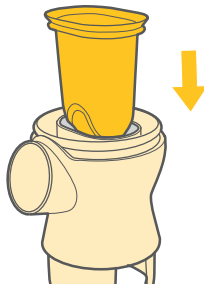2

- Carefully insert the (dark yellow) membrane into the top of the connector body.
- The membrane must be secure for correct assembly.

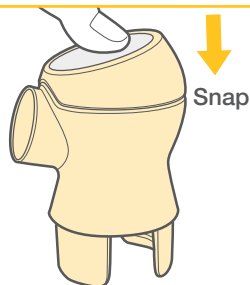3


Press the connector cap onto the connector body. It will snap closed.

5


Fasten the bottle onto the complete breast shield assembly.

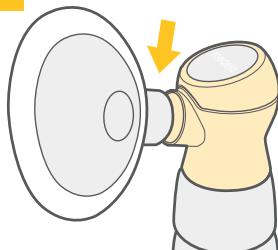6


Push breast shield onto the connector body.

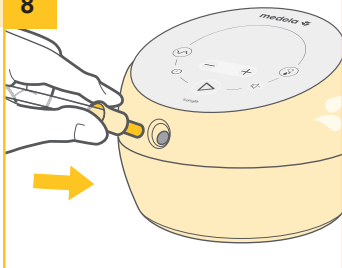8

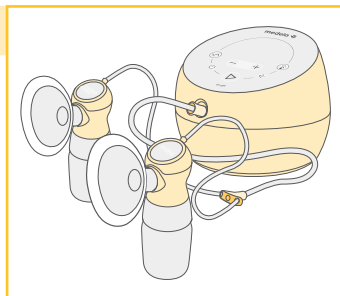Insert the long tubing end into the breast pump as far as it will go.


Accurately assembled system.

10. Powering Your Breast Pump

10A. Charging the Breast Pump

There are two ways to power your breast pump for your pumping session:

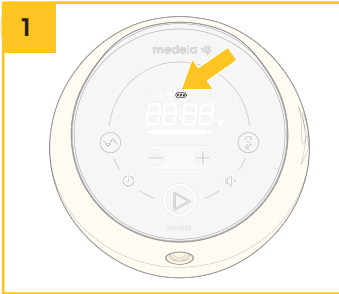
Rechargeable Battery

- Before first use, charge the breast pump battery for 12 hours with the power adaptor. The battery icon will appear fully charged on the display when charging is complete.
- Proceed to **Section 11**, Operating Your Breast Pump, to begin your pumping session.

OR

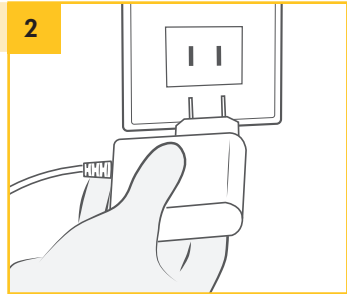Power Adaptor

- Plug power cord into the power adaptor port on the back of the breast pump. Plug power cord into power source.
- Proceed to **Section 11**, Operating Your Breast Pump to begin your pumping session.

10B. Battery Care & Status Codes

Your breast pump contains a Lithium ION rechargeable battery.

T Tip

- Before using on battery power for the first time, fully charge your breast pump for 12 hours.
- You can use your breast pump while charging the battery. The battery will charge as soon as the breast pump is connected to a power source. It will not charge while the pump is running.
- Keep your pump in a cool place.
- Keep your battery charged.
- Fully charge the battery before storing your pump for an extended period of time.
- Recharge the battery before it completely drains. This is better for the lifecycle of the battery.
- If your breast pump has been stored in a hot location, it may not run on battery power right away. To resume normal battery function, allow the breast pump to cool for one hour. During this time, you can use your breast pump when plugged into external wall power.


Battery Charge Status Codes

	External power source connected
	Battery charge status
 (flashing)	Low battery, needs charging
 bAtt	Battery drained. See Troubleshooting (Section 16).
 (flashing)	Battery broken. See Troubleshooting (Section 16) and Call Medela Customer Service.



10. Powering Your Breast Pump (cont.)

10C. Power Adaptor Care

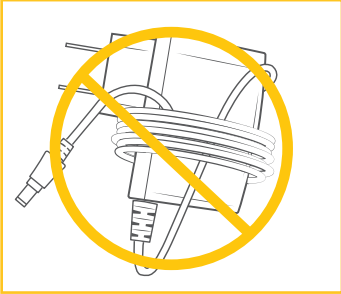
How to care for your adaptor:

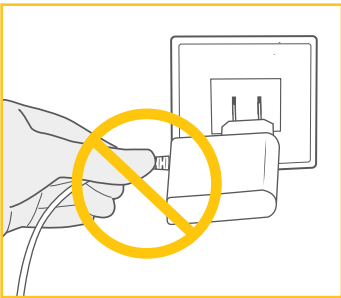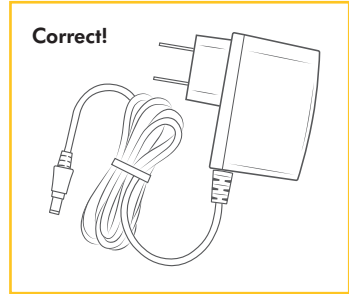Do not wrap cord of power adaptor around the plug body.

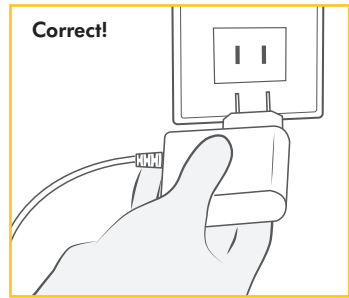Do not unplug power adaptor by pulling on the cord.

Never operate an electrical device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. If damage is found, immediately discontinue use of power adaptor and call Medela Customer Service at 1-800-435-8316.

Warning

- Use only the power adaptor that comes with the breast pump.
- Make sure the voltage of the power adaptor is compatible with the power source.

10D. Traveling Outside of Your Country



When travelling internationally, we recommend purchasing a Universal Power Plug Adaptor (not included). Please consult with the country you are visiting to find out what type of adaptor will work best. This breast pump's A/C adaptor is two-pronged and does not require a grounded outlet. This breast pump has a Lithium ION Battery which may be restricted from the country you are visiting. Please consult with the country you are visiting to find out if there are any restrictions that pertain to travelling with Lithium ION.

For information about traveling with your breast pump and expressed breast milk, visit www.TSA.gov.

11. Operating Your Breast Pump


11A. Glossary of Terms, Controls & Icons

<p>2-Phase Expression®</p> 	<p>Through research, Medela found that babies nurse in 2 Phases – Stimulation and Expression. This research is the basis for the technology in all Medela breast pumps.</p>
<p>Stimulation Phase</p>	<p>Fast sucking/pumping rhythm to stimulate milk flow.</p>
<p>Expression Phase</p>	<p>Slower sucking/pumping rhythm for gentle and efficient milk removal after milk has started flowing.</p>
<p>Maximum Comfort Vacuum™</p>	<p>Highest vacuum setting where pumping still feels comfortable. Different for every mother.</p>
<p>Pumping Rhythms</p>	<p>This breast pump includes two clinically-tested pumping rhythms. It offers you options to choose the most comfortable and effective for you.</p>
<p>Responsive Pump Technology</p>	<p>Responsive Pump Technology automatically adjusts to the body and environment and makes sessions more consistent and productive.</p>
<p>Letdown</p>	<p>When milk begins to flow from the breast.</p>

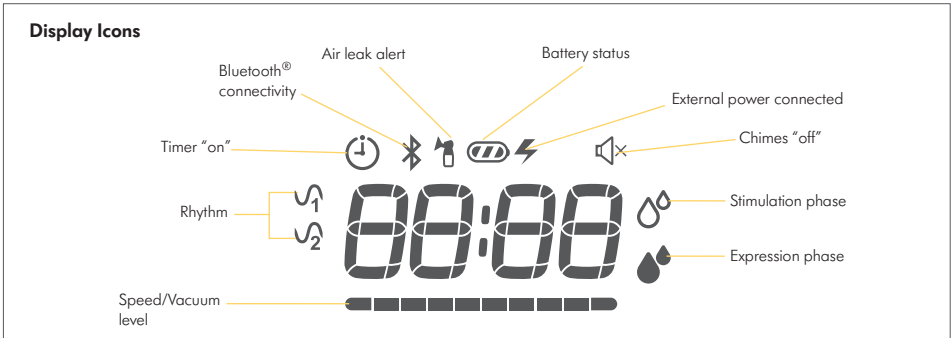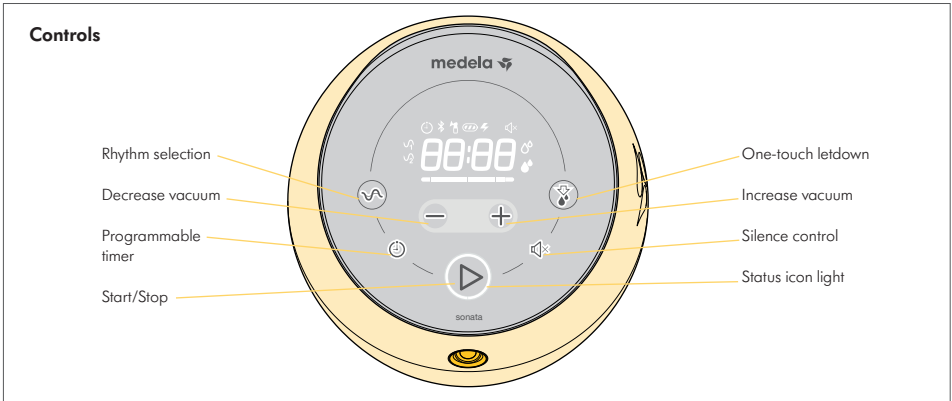
11. Operating Your Breast Pump (cont.)

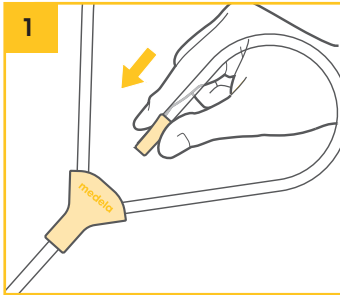11B. Single Pumping

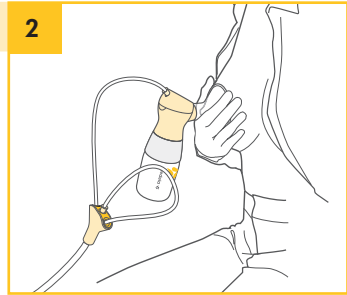
1
Insert the unused tubing end into the tubing holder.

2
Hold the breast shield and connector onto your breast with your thumb and index finger. Support your breast with the palm of your hand.

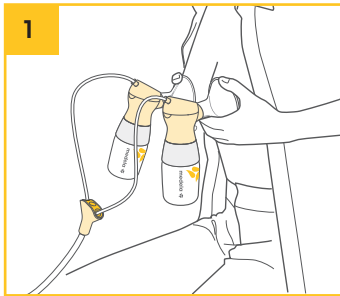11C. Double Pumping


1
Place the breast shields on your breasts so that your nipples are properly centered in the tunnels.

Refer to Section 8 for breast shield sizing information.


T Tip

- Wash hands thoroughly with soap and water before touching breast pump, kit and breasts. Avoid touching the inside of containers or lids.
- Always inspect all parts prior to use for cleanliness.

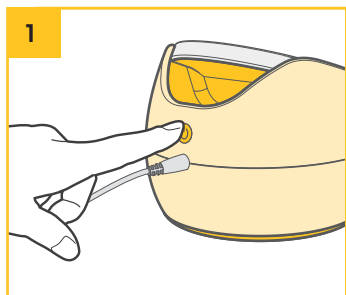
11D. Pumping Steps

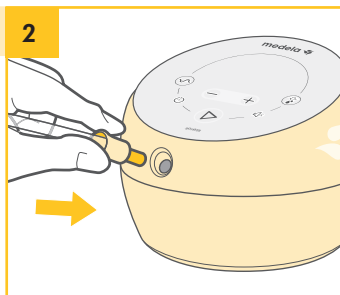Turn on the breast pump by pressing the power button  on the back.

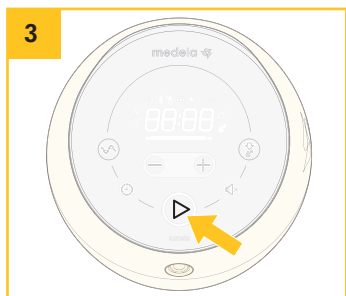Plug your assembled double-breast pump kit into the tubing port on the front of the breast pump.


Start breast pump vacuum by pressing the start/stop  button.

Caution


- Make sure tubing is not kinked or pinched while pumping.
- Do not hold the breast pump kit by the bottle. This can lead to blockage of the milk ducts and engorgement.
- Do not try to express with vacuum that is too high and uncomfortable (painful). The pain, along with potential breast and nipple trauma, may decrease milk output. The comfortable level is based on each individual.
- Contact your healthcare professional or breastfeeding specialist if you can express only minimal or no milk or if expression is painful.

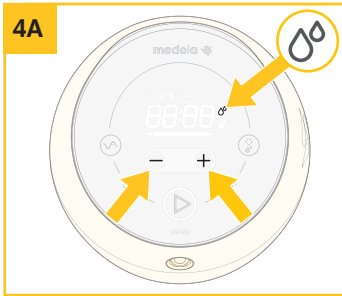11. Operating Your Breast Pump (cont.)

Your breast pump will begin in Stimulation Phase. Adjust the pumping level using the controls to find a comfortable level.

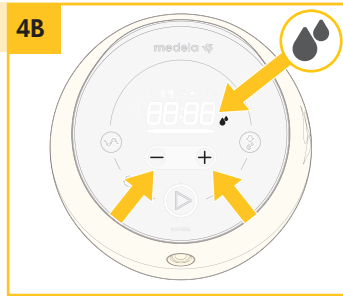
Stimulation Phase – a fast sucking/pumping rhythm to stimulate milk flow.


After two minutes, your pump will automatically transition to Expression Phase. Adjust the pumping level using the controls to find your Maximum Comfort Vacuum **(see Section 12).**

If your milk begins to flow before two minutes, press the let-down button manually transition to Expression Phase.

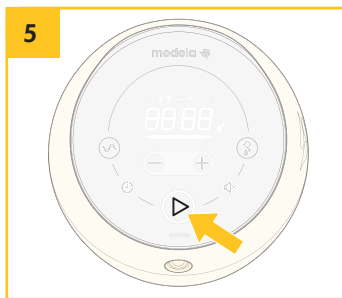Expression Phase - a slower sucking/pumping rhythm for gentle and efficient milk removal after milk has started flowing.

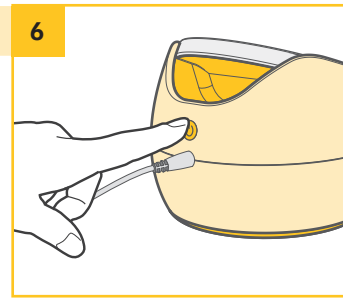When your pumping session is over, stop vacuum by pressing the start/stop button .

To turn off your breast pump, press the power button on the back.

T Tip

- The breast pump will automatically turn off after five minutes of inactivity.
- Do not forget to record milk output in your MyMedela app.

12. Finding Your Maximum Comfort Vacuum™

Once you are pumping in the **Expression Phase**,  increase the speed/vacuum \oplus until pumping feels slightly uncomfortable (not painful), then decrease \ominus slightly.

Tip

- Maximum Comfort Vacuum is the highest vacuum setting where pumping still feels comfortable. This is different for every mother.
- Stimulation should be at a comfortable vacuum level.
- Reassess your Maximum Comfort Vacuum throughout your pumping experience. It can change throughout each stage of lactation.

13. Features

Breast Pump Connectivity

By tracking your sessions in our MyMedela® app, you can stay on-target with daily highlights and custom content. Having visibility to your pumping history gives you the opportunity to plan efficiently each day and track changes in your body like an increase or decrease in milk output.

By connecting your breast pump to your mobile device you can automatically transfer your pumping data (session length, phases & levels) and manually input the amount of milk expressed into your MyMedela pumping log.

How do I connect (pair) my Sonata Breast Pump with my mobile device?

Step 1. Download and install MyMedela  to your mobile device. Available for free for iOS and Android.

Step 2. Open MyMedela on your mobile device. Follow set-up and pairing instructions.

Step 3. When correctly paired, the  will appear on your breast pump display.

**Tip**

- After initial set-up, the Bluetooth® icon  will automatically appear on your breast pump, when powered on, if your paired device is near and MyMedela is open.
- You can pair multiple mobile devices with Sonata.
- Only one mobile device at a time can interact with Sonata.
- Sonata has the ability to store your last 30 pump sessions in its memory. If the pump is not connected prior to session 31 your history will be overwritten with new session data. When connected, MyMedela will log an infinite number of historical sessions.

How do I remove my Sonata breast pump from my mobile device?

Step 1. Navigate through MyMedela to remove your Sonata Breast Pump.

**Tip**

- Removing Sonata from your mobile device will not delete information stored in MyMedela.
- If you re-connect your breast pump and mobile device the last 30 records stored within your pump will transfer.

To learn more about MyMedela, including download instructions, please visit:

www.mymedelaapp.com.

13. Features (cont.)

30

Selecting a Rhythm

Two clinically tested pumping rhythms are included:

1. **Signature** - This rhythm is based on Medela's successful Symphony® Breast Pump used in hospitals.
2. **Lifestyle** - An alternate rhythm for when you want to change your pumping routine. It includes additional benefits of longer battery, and quieter pump operation when you need it.

Note

Your breast pump will default to Signature rhythm. To change to the alternate Lifestyle rhythm, push the pattern  button.

Setting the Session Timer

The default timer on your breast pump will start at zero and count up. Follow the steps below if you want to count down from a set time.

To set your count-down timer:

1. Before you begin pumping, press the timer  button. The timer icon will appear on the display and the   buttons will be illuminated.
2. Use the   buttons to adjust time in 1 minute increments (up to 30 minutes).
3. Press start  to begin your pumping session.
4. Once you've reached your goal time, the vacuum will pause and the display will show your total pumping duration.

Note

- Your breast pump will make a sound and the display will flash to let you know that your session has 1 minute remaining.
- Your breast pump will make a sound and flash to alert you that your session has ended.

Tip

If you want to continue the same session, press start/stop  once and vacuum will continue in same level, phase and rhythm. The timer will count up from your previously entered pumping time.

13. Features (cont.)

32

Pausing Your Breast Pump

Your breast pump has the ability to briefly pause during your pumping session. This allows you to readjust yourself or your breast pump pieces, or attend to your surroundings.

To pause your breast pump:

1. Press the start/stop  button. Vacuum will stop and time on the display will flash. Start/stop  icon on display will flash.
2. When ready to resume your session, press start/stop  button.

Note

- If pumping is not resumed within two minutes, the breast pump will end your session. To resume pumping, you will need to start a new session.
- If you need to pause for more than two minutes, we recommend restarting in Stimulation Phase.

Silencing Your Breast Pump Chimes

Your breast pump is equipped to provide sound notifications, such as button presses, during your pumping sessions.

To silence your breast pump chimes:

1. Push the Silence control  button.
2. Your breast pump will remember your setting for future sessions. To reactivate the sound feature, push the button again.

Note

When the breast pump chimes are turned off, the sound  icon will appear on the display screen.

14. Storing Breast Milk

Storing Breast Milk in a Medela Cooler Bag

- Your Medela cooler bag requires one Medela contoured ice pack. Up to four 5 oz/150 mL bottles of breast milk will stay cool for up to 10 hours in a room temperature environment inside your pump bag. Up to four 5 oz/150 mL bottles of breast milk will stay cool for up to 8 hours in a room temperature environment, inside the cooler, stored outside of the pump bag.
- As soon as you are finished pumping, place the breast milk bottle or breast milk storage bag in the cooler bag with the ice pack.
- Transfer breast milk bottles or breast milk storage bags to the refrigerator or freezer once you are home.

Storage Guidelines for Freshly Expressed Breast Milk (For Healthy Term Babies)			
Room Temperature 60°F to 77°F (16° C to 25° C)	Refrigerator 39°F (4°C) or colder	Freezer 0°F (-18°C) or colder	Breast Milk Thawed in the Refrigerator
Up to 4 hours is best *Up to 6 hours for milk expressed under very clean conditions	Up to 3 days is best *Up to 5 days for milk expressed under very clean conditions	Up to 6 months is best *Up to 9 months for milk expressed under very clean conditions	At room temperature: Up to 2 hours Refrigerator: Up to 24 hours DO NOT REFREEZE

References: www.BreastmilkGuidelines.com

15. Preparing & Feeding Breast Milk

Preparing Breast Milk

Caution

- Do not thaw frozen breast milk in a microwave or in a pan of boiling water.
 - Do not microwave breast milk. Microwaving can cause severe burns to baby's mouth from hot spots that develop in the milk during microwaving. (Microwaving can also change the composition of breast milk.)
-
- Thaw breast milk overnight in the refrigerator. Thawed breast milk is safe in the refrigerator for 24 hours. Do not refreeze thawed breast milk.
 - Quickly thaw breast milk by holding the bottle or bag under warm running water.
 - Place the sealed bottle or bag in a bowl of warm water for 20-minutes to bring it to body temperature.
 - If you are adding expressed breast milk to a container of already frozen breast milk, make sure to add a lesser amount than the already frozen amount.

Feeding Breast Milk

It is recommended that breastfeeding is well established prior to bottle feeding your baby.

- Always inspect the bottle, nipple and other components immediately before and after each use. If nipple appears cracked or torn, discontinue use immediately.
- To prevent possible choking hazard, test strength of nipple by pulling on bulb portion of the nipple.
- Do not attempt to enlarge the nipple hole.
- Infants must not be bottle fed without adult supervision.
- Nipple should not be used as a pacifier.

For additional breast milk collection & storage information, please visit www.BreastmilkGuidelines.com.

16. Troubleshooting

Air Leak Alert

There is an air leak in your system and vacuum is not operating at optimal performance. During this time, the breast pump will continue to run, but an alert icon will appear on the display and a chime will sound every 30 seconds until the correction has been made.

- Verify that the kit connections are properly assembled.
- Make sure the breast shields are secure to your breasts.
- Check that the tubing connection between the breast pump and kit is secure.
- If you still have issues after following these steps, please call Medela Customer Service at 1-800-435-8316 for further assistance.

Note

There is no need to turn off your breast pump. The alert  icon will disappear once your breast pump detects that the error has been fixed. In order to do this, the breast pump must be running.

Battery Drained

If bAtt appears on the screen with an empty battery cell, the battery is too drained to start a vacuum.

To recover:

- Plug your breast pump into an external power source.
- Fully charge the battery for 12 hours.

16. Troubleshooting (cont.)

38

Battery Broken

If the battery icon is empty and flashing, your battery is broken. The battery cannot be serviced or replaced. Please call Customer Service at 1-800-435-8316 for replacement options.

To continue pumping:

Plug your breast pump into an external power source to use.

Note

If your battery is broken, your breast pump will no longer be able to charge or run on battery power.

Bluetooth® Broken

If Err flashes on your breast pump for three seconds it means the Bluetooth® connection is broken and the breast pump or mobile device are unable to complete data transmission. This will not impact the ability to pump with Sonata, only the ability to transfer data to the MyMedela® app. Please call Customer Service at 1-800-435-8316 for assistance.

16. Troubleshooting (cont.)

40

General Alert

Pump is not operable in this condition. Please contact Medela Customer Service at 1-800-435-8316 for assistance.

Additional Troubleshooting

<p>If your breast pump is not responding as expected</p>	<ul style="list-style-type: none"> To reset your device, plug Sonata into external power and press and hold the power button (on the back) for up to 10 seconds or until the pump display flashes. If there is no change contact Medela Customer Service at 1-800-435-8316.
<p>If breast pump does not power on</p>	<ul style="list-style-type: none"> Check that the breast pump is plugged into a working electric outlet. If there is no power, display does not light up, and there is no error alert appearing on screen, contact Medela Customer Service at 1-800-435-8316.
<p>If breast pump does not run after pushing the start/stop button </p>	<ul style="list-style-type: none"> Make sure you've powered on the breast pump. Make sure breast pump is plugged in and/ or charged. For more information, refer to General Alert (Section 16).
<p>If breast pump exterior gets wet</p>	<ul style="list-style-type: none"> Unplug from power source and power off and discontinue use. Dry off the exterior of the breast pump.
<p>If breast pump is submerged in water</p>	<ul style="list-style-type: none"> Unplug from power source, power off and discontinue use. Contact Medela Customer Service at 1-800-435-8316.
<p>If there is low or no suction</p>	<ul style="list-style-type: none"> Check that all of the breast pump and kit connections are secure. Disassemble the double-pumping kit and check that all pieces are clean and not damaged. Assemble the pump set as directed in Assembling your Breast Pump (Section 9). While pumping, make sure the breast shields form a complete seal around the breast. When single pumping, check that the unused end is correctly stored in the tubing holder. If suction does not improve after following these steps, contact Medela Customer Service at 1-800-435-8316.
<p>If there is a breast milk overflow</p>	<ul style="list-style-type: none"> Turn off the breast pump and unplug the power cord from the electrical outlet. Use a damp (not wet) cloth to wipe the breast pump. If there is milk in the tubing: Remove tubing from the pump unit as well as breast shield assembly. Clean as directed in Section 6. When all components are dry, reassemble.
<p>No "let-down" or breast milk expression</p>	<ul style="list-style-type: none"> Ensure your breast pump kit is assembled correctly and suction is available. Relax and take a 10-15 minute break if let-down is not achieved after two consecutive pumping sessions. Consult with a health care or lactation professional if expression does not occur.
<p>Data cannot be sent</p>	<ul style="list-style-type: none"> The Bluetooth® function  on the mobile device is turned off. Turn on the Bluetooth® function and try sending data again. Pairing between the two devices has not been completed. Perform pairing instructions again. Refer to "How do I pair my Sonata Breast Pump with my mobile device?" in Section 13. The mobile device is too far away from the breast pump. Move the device and breast pump closer together. The application on the mobile device is not ready. Check the application. If the Err symbol appears after checking the application, contact Customer Service at 1-800-435-8316.

If you have not resolved the problem with your breast pump or you have further questions, please contact Medela Customer Service at 1-800-435-8316 or fill out the Contact Us form at www.medelabreastfeedingus.com | www.medela.ca.

Questions? Visit www.medela.com or call us at 1-800-435-8316.

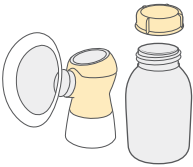
17. Long-term Storage

42

Your breast pump contains a rechargeable battery. To preserve your breast pump, when not in use for an extended amount of time, it will automatically transition into storage mode.

To remove from storage mode, plug your breast pump into an external power source and fully charge the battery for up to 12 hours.

✓ **Before storing your breast pump for an extended amount of time remember to:**

Clean breast pump kit
(Section 5)

Clean breast pump
(Section 6)

Make sure breast pump battery is fully-charged

18. Disposal

Disposal

The symbol on the product or its packaging indicates that this product must not be disposed of with your other household waste. Instead, it is your responsibility to dispose of your waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste for recycling, please contact your local authority, or where you purchased your product.

Recycle your breast pump

Through the Medela Recycles program, Medela LLC. provides customers the option to recycle their breast pump after they have finished their breast milk feeding journey. For information on how to recycle your pump, visit www.MedelaRecycles.com.

This program is currently only available to U.S. residents. We only accept the return of electric breast pumps. All other breast pump kit components (connectors, breast shields, tubing, bottles, valves, membranes, and cooling elements) can be disposed of in your waste management recycling container at your home. Cooler bags and other soft-good items are not recyclable.

19. Supplemental Information

Supplemental Information

Following are some common breastfeeding-related conditions. If you experience any of these symptoms, contact a healthcare professional or lactation professional.

	Symptom	Potential Causes
Engorgement	Breasts are hard and uncomfortable, possibly with reddened areas.	Milk not draining properly from breast. Milk can drain back into the tissue causing areas to swell and become tender.
Blocked/Plugged Ducts	Area of breast will look reddened and may be tender to touch.	Milk not being drained from a specific duct. The area becomes "clogged" up and milk is then prevented from flowing.
Mastitis	Often follows engorgement but can just occur suddenly. The initial symptoms are similar to flu – fatigue, headache and muscle aches, fever and localized breast tenderness. It can be in one or both breasts and needs immediate treatment.	A bacterial infection in the breast tissue which is often accompanied by cracked nipples.

A Healthcare Professional or Lactation Professional can give you guidance in breastfeeding your baby and using a breast pump. Their guidance is valuable in terms of the lifetime of health benefits for both you and your baby. A listing of lactation professionals are available at www.medelaLocator.com or calling 1-800-TELL-YOU (1-800-835-5968).

20. Warranty

This product is warranted by Medela to the original retail purchaser to be free from defects in material and workmanship for the period of 1 year for pump mechanism (90 days for other parts and accessories) from the date of purchase. Warranty can only be claimed in the country of purchase. In the event of defect, at Medela's option, Medela will replace this product, without charge for such replacement parts. Purchaser shall bear all expense for returning this product to Medela. This warranty does not apply to any product used commercially or which has been subjected to misuse, abuse or alteration.

ANY AND ALL IMPLIED WARRANTIES, INCLUDING THE WARRANTY OF MERCHANTABILITY, ARE LIMITED TO A DURATION OF 3 YEARS FROM DATE OF PURCHASE. SOME AREAS DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE LIMITATIONS MAY NOT APPLY TO YOU. THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS AND YOU MAY ALSO HAVE OTHER RIGHTS WHICH VARY FROM STATE TO STATE OR IN YOUR COUNTRY.

Before you make a claim under this warranty, it may save time and expense to call Medela Customer Service (toll free) at 1-800-435-8316. You may also call this number for additional information concerning this warranty. If you wish to make a claim under this warranty, you must return this product to Medela with a return authorization number received from Medela Customer Service, prepaid, together with your dated bill of sale or other proof of purchase and a brief statement of the problem to the following address:

Medela LLC – Returns, Door 4501
1101 Corporate Dr.
McHenry, IL 60050
ATTENTION: RETURNS

Medela Canada
4160 Sladeview Crescent Unit #8
Mississauga, ON L5L 0A1
ATTENTION: RETURNS

Contact Customer Service for a return authorization number. Returns are not accepted without an authorization number.

We recommend recording your serial number here when you first open your pump for easy reference when contacting Medela Customer Service.

Your serial number can be found on the sticker on the bottom of the breast pump.

SERIAL NUMBER:

21. Meaning of Symbols

Glossary of Symbols

 On/off

 Session duration

 Refer to Operating Instructions for important safety information, failure to follow instructions could place operator at risk. (when shown on the product this symbol is blue)

 This symbol indicates the manufacturer.

 This symbol indicates do not dispose the device together with unsorted municipal waste (in accordance with local regulations).

 This symbol indicates compliance with international requirements for protection from electric shock. (Type BF applied parts.)

IP21 This symbol indicates the protection against entrance of solid foreign objects and against harmful effects due to the entrance of water.

 This symbol indicates the date of manufacture (four digits for the year and two digits for the month).

 This symbol indicates Class II.

 This symbol indicates that the power adaptor is for indoor use only.

 Classified by UL according to the standards ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014), IEC 60601-1-6:2010 + AMD1(2013), IEC 62366:2007 + AMD1(2014) and IEC 60601-1-11:2015

 This symbol indicates alternating current.

 This symbol indicates direct current.

CE Symbol showing that the manufacturer has taken all necessary measures to ensure that the product complies with the applicable safety legislation for distribution within the European Community. This symbol only applies to the power adaptor.

 This device contains Bluetooth® wireless technology. “The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc., and any use of such marks by Medela LLC. is under license. Other trademarks and trade names are those of their respective owners”.

 Warning sign (when shown on product this symbol is yellow)

 UL Recognized Component.

 Keep Dry.

FC Certification mark on electronic products manufactured or sold in the United States which certifies that the electromagnetic interference from the device is under limits approved by the Federal Communications Commission.

22. EMC/Technical Description

The breast pump needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the instructions for use. Portable and mobile RF communications can affect the breast pump.

Note

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

22. EMC/Technical Description (cont.)

Guidance and manufacturer’s declaration – electromagnetic emissions

This breast pump is intended for use in the electromagnetic environment specified below. The customer or the user of the breast pump should assure that it is used in such an environment.

Emission tests	Compliance	Electromagnetic environment – guidance
RF Emissions CISPR 11	Group 1	The breast pump uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions	Class B	The breast pump is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Pst < 1.0	

Warning

Warning – This breast pump should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, this breast pump should be observed to verify normal operation in the configuration in which it will be used.

Guidance and manufacturer's declaration – electromagnetic immunity

This breast pump is intended for use in the electromagnetic environment specified below. The customer or the user of the breast pump should assure that it is used in such an environment. The Sonata Breast Pump has no essential performance but was tested for immunity to electromagnetic disturbances and passed using the following criteria:

1. No visible change in the operation of the breast pump.
2. The breast pump changes settings but returns automatically to previous settings.
3. The breast pump changes settings but can return to previous settings by intervention of the user.
4. The breast pump changes settings but can return to previous settings by executing the INITIATION SEQUENCE.

INITIATION SEQUENCE: (The following sequence can be used to reset the device in the event of an EM Disturbance)

1. Unplug the power cord from the back of the breast pump.
2. Unplug the power adaptor from the wall outlet.
3. Attach a breast pump kit to the pump port (if none is attached).
4. Press and hold the power button for at least 30 seconds. (Note: The breast pump may cycle on/off several times. After 30 seconds the breast pump should remain off.)
5. Plug the power adaptor back into the wall outlet.
6. Plug the power cord into the back of the breast pump.
7. Press the on/off switch to power on the breast pump.
8. Press the 'PLAY' button to begin pumping.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	+/- 2kV, +/- 4kV, +/- 6 kV, +/- 8kV contact discharge +/- 2kV, +/- 4kV, +/- 6 kV, +/- 8kV, +/- 15 kV air discharge	+/- 2kV, +/- 4kV, +/- 6 kV, +/- 8kV contact discharge +/- 2kV, +/- 4kV, +/- 6 kV, +/- 8kV, +/- 15 kV air discharge	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV 100 kHz repetition frequency	± 2 kV 100 kHz repetition frequency	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0,5 kV, ± 1 kV Line-to-line	± 0,5 kV, ± 1 kV Line-to-line	Mains power quality should be that of a typical commercial or hospital environment.

22. EMC/Technical Description (cont.)

<p>Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11</p>	<p>0 % UT ; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°</p>	<p>0 % UT ; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°</p>	<p>Mains power quality should be that of a typical commercial or hospital environment. If the user of the breast pump requires continued operation during power mains interruptions, it is recommended that the breast pump be powered from an uninterruptible power supply or a battery.</p>
<p>Power frequency (50/60 Hz) magnetic field IEC 61000-4-8</p>	<p>30 A/m, 50 or 60 Hz</p>	<p>30 A/m, 50/60 Hz</p>	<p>It may be necessary to position the breast pump further from sources of power frequency magnetic fields or to install magnetic shielding. The power frequency magnetic field should be measured in the intended installation location to assure that it is sufficiently low.</p>

NOTE U_T is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration – electromagnetic immunity

This breast pump is intended for use in the electromagnetic environment specified below. The customer or the user of the breast pump should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz	6V	<p>Portable and mobile RF communications equipment should be used no closer to any part of the breast pump, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommend separation distance $d = 1,2\sqrt{P}$ $d = 1,2\sqrt{P}$ 80 MHz - 800 MHz $d = 2,3\sqrt{P}$ 800 MHz - 2.7 GHz</p> <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol.</p> 
Radiated RF IEC 61000-4-3	10 V/m (minimum) 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m (minimum)	

22. EMC/Technical Description (cont.)

Note 1 At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the breast pump is used exceeds the applicable RF compliance level above, the breast pump should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the breast pump.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the breast pump

This electric breast pump is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the electric breast pump can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the electric breast pump as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800 MHz to 2,7 GHz $d = 2.3\sqrt{P}$
0,01	0,06	0,04	0,07
0,1	0,18	0,11	0,22
1	0,58	0,35	0,22
10	1,84	1,11	2,21
100	5,83	3,50	7,00

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the breast pump including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

23. Technical Specifications

Technical Specifications

Vacuum Range
-50.....250 mmHg
54 ... 120 cpm

Power In/Power Out



In	Out
100-240V- 50/60 Hz	12V --- 1.5 A ---

Battery Capacity and Type
7.4 V, 2200mAh
Li-Ion

Size
160 x 160 x 120

Weight
1.1 kg

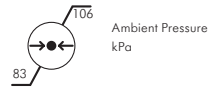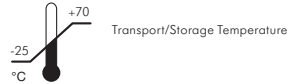Bluetooth® RF Output Power 0dBm (class 3)
Frequency band 2379–2496 MHz Frequency-hopping spread spectrum per Bluetooth® 4.0 Core Specification.

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info@medela.ca

www.medela.ca

Content appearance may vary from pictures.

Medela wordmark and logo, Sonata, Symphony, 2-Phase Expression, MyMedela and Easy Expression are registered in the U.S. Patent and Trademark Office. Mother's Milk, Everyday Amazing, Responsive Pump Technology, Quick Clean, Micro-Steam, Maximum Comfort Vacuum, PersonalFit and PersonalFit Flex are trademarks of Medela.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Medela is under license.

iOS is a trademark of Cisco in the U.S. and other countries. Medela is not affiliated with Cisco. Android is a trademark of Google Inc. Medela is not affiliated with Google.

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Medela Canada, Inc.